YE-690A

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yuwell

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
1. PRODUCT FEATURES

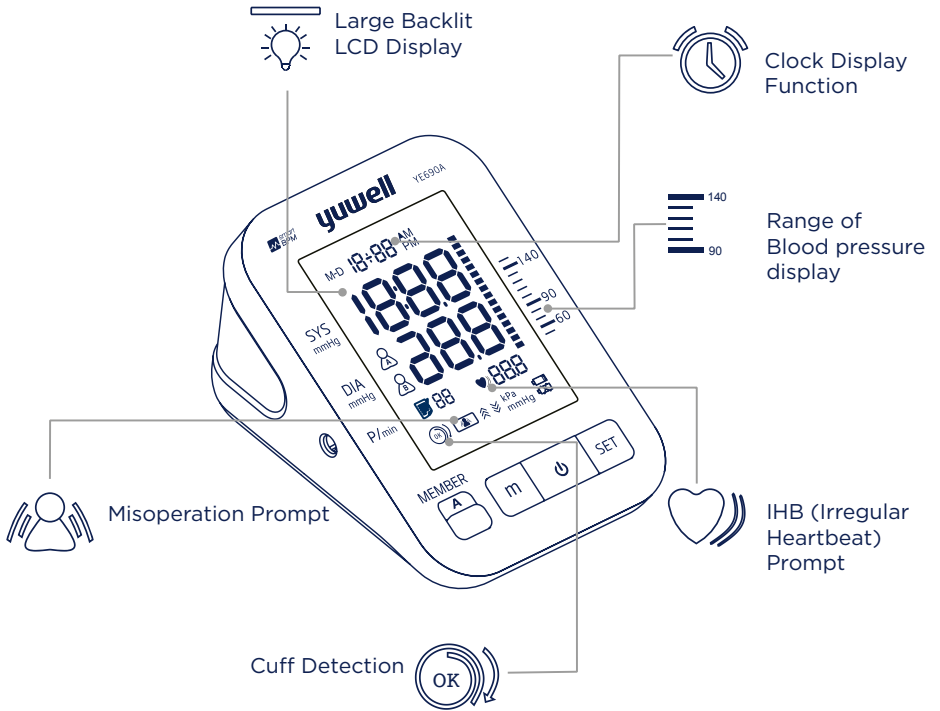
2. OPERATION GUIDE


Please follow the steps as shown below

2.1. POWER-ON


Connect the power source

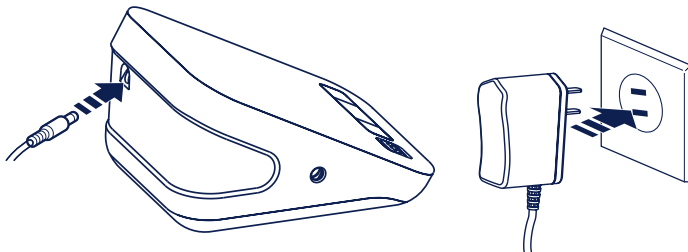
Or install four AA batteries

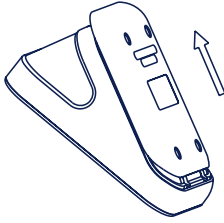Open the battery cover

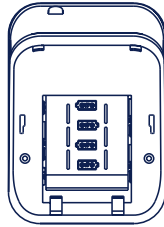Install the batteries properly


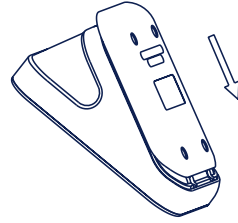
Close the battery cover

(Please turn to page 9 for details)

2.2. MEASUREMENT


1 USING METHOD OF CUFF:
Wearing the cuff and make the lower edge 2-3cm above the el-

bow, keep the air tube inside of the arm.

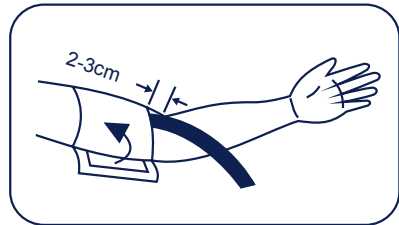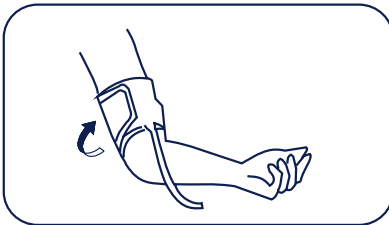2 MEASURING POSTURE
Sit straightly, keep the center of

cuff and heart at the same level.

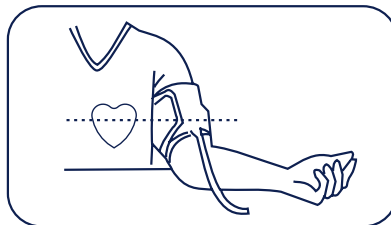3 START MEASURING
Sit still for 5 minutes and then Press the „  „ Button to

measure (Please turn to page 18 for details)


2.3. CHECK THE RECORDS

Press the button of „ **m** „, to check the measuring records.
This monitor can display the ave-

rage measurement and store 74× 2 records. (Please turn to page 29 for details)


3. PRODUCT STRUCTURE

1 Main part

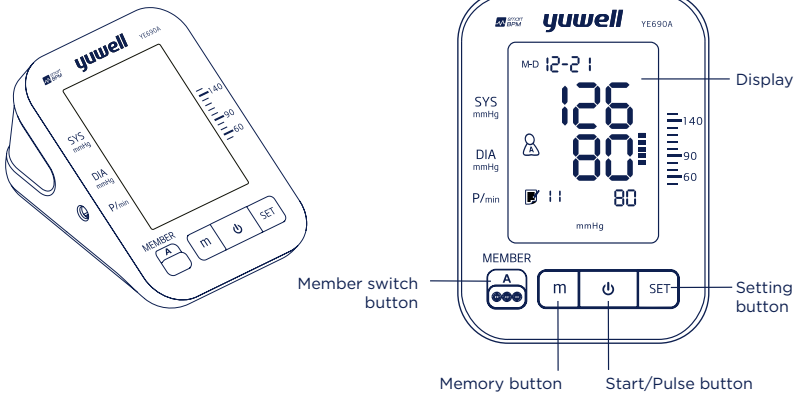
2 Mankiet


The suitable upper-arm circumference of the cuff is 22cm-45cm

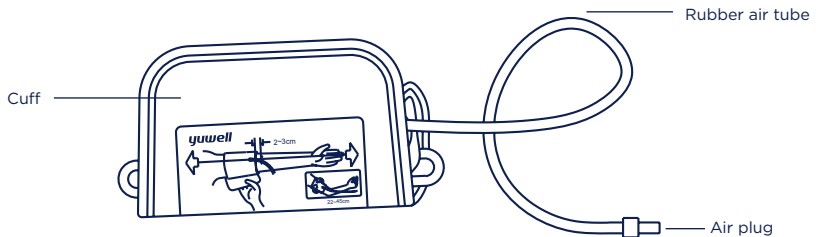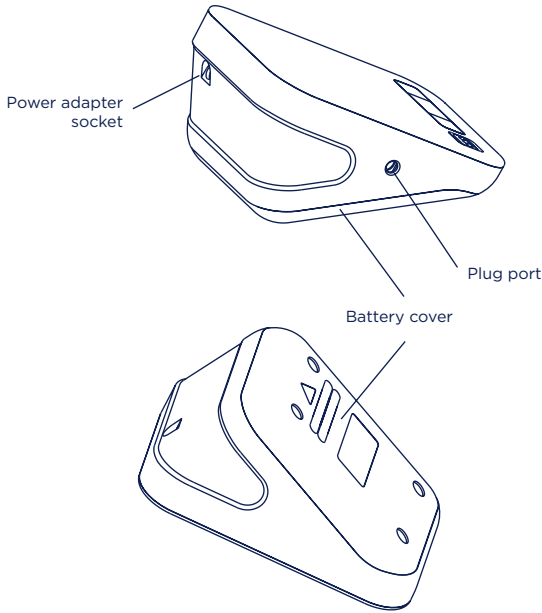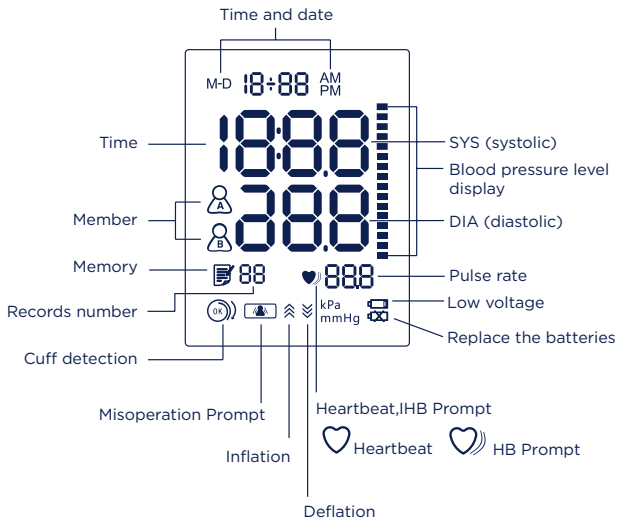3 Power

4 Display

3. SYMBOLE

SYMBOL	DESCRIPTION	SYMBOL	DESCRIPTION
	Attention!!		Keep Upright
	Class II device		Fragile
	Type BF application part		Keep Dry

4. APPLICATION SCOPE AND ATTENTION ITEMS

4.1. RANGE OF APPLICATION

The monitor is used for measuring blood pressure of human body (adu-

It older than 12, not suitable for new born baby)!

4.1. WSKAZÓWKI I ZALECENIA

CAUTION: PAY ATTENTION TO THE FOLLOWING POINTS WHEN MEASURING OR IT MAY CAUSE THE INCORRECT RESULTS.

SIT STILL FOR 5 MINUTES BEFORE MEASURING TO ENSURE QUIET AND STABLE MODE.

DO NOT TAKE THE MEASUREMENT WHILE STANDING, WALKING OR HAVING BODY PRESSED.

DO NOT TAKE THE MEASUREMENT AFTER SMOKING, DRINKING WINE OR COFFEE (BLACK TEA).

DO NOT TAKE THE MEASUREMENT AFTER SPORT OR BATH.

DO NOT SPEAK OR MOVE WHILE MEASURING.

DO NOT MOVE, SHAKE ARM OR BEND FINGERS WHILE MEASURING.

DO NOT TAKE THE MEASUREMENT AT EXTREME TEMPERATURE CONDITION OR THE VARIED SEVERELY ENVIRONMENT.

DO NOT TAKE THE MEASUREMENT IN A MOVING VEHICLE.

DO NOT MEASURE CONTINUOUSLY. (2-3 MINUTES OR MORE SHOULD BE SPARED BETWEEN TWO MEASUREMENTS).


DO NOT MEASURE WITHIN 1 HOUR AFTER MEAL.
 DO NOT USE MOBILE PHONE NEAR THE MONITOR.
 PLEASE REINSTALL THE BATTERIES AND START AGAIN IF CAN-
 NOT MEASURE.
 FOR PATIENT OF ARRHYTHMIA, MEASURING RESULTS MAY NOT
 BE ACCURATE.
 DO NOT KEEP THE CUFF IN THE AERATED STATE FOR A LONG
 TIME.
NOTE: DO NOT DIAGNOSE WITH THE MEASUREMENT, PLEASE
 FOLLOW DOCTOR'S INSTRUCTION.


STATEMENT: IF THE MONITOR HAS NOT BEEN STORED IN THE
 REQUIRED TEMPERATURE AND HUMIDITY RANGE, IT MAY NOT
 CONFORMANCE TO SPECIFICATION.

4. BATTERY INSTALLATION


- Open the battery cover and install four AAA batteries properly.
- Then close the cover with a click sound.

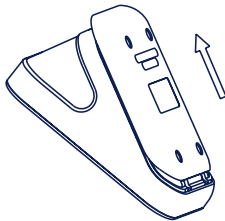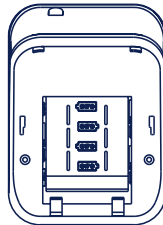Open the battery cover

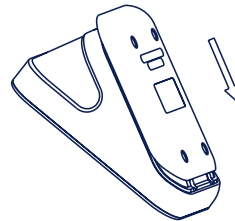Install the batteries properly

Close the battery cover

- The "  „ icon appearing means the battery is low.
- The "  „ icon appearing means the battery is running out. Please replace all the batteries.
- Please take out the batteries if the monitor will not be used for a long time (over three months).
- Please take out the batteries if use the DC external power for a long time.

6V/600MA DC EXTERNAL POWER CAN BE CONNECTED TO THE
 MONITOR. (PLEASE USE THE POWER ADAPTER SUPPLIED BY US.)

5. TIME AND DATE SETTING

The monitor will store the measuring results automatically by the management system.


It is necessary to reset time and date after installing new batteries or connecting to DC power. Please operate as following steps. (For example: setting the date as 2015-1-15 and time as 08:28)

1. Year setting: Press the "  ", button for more than 3 seconds till the number starts flashing.
2. The year increase once press the "  ", button, and decrease once press the "  ", button.
3. Press the "  ", month setting.
4. Using the same way for other settings.

6. MMHG/KPA SETTING

After finishing the date and time setting, pressing the "  ", button to enter the mmHg/ kPa unit setting. Press the "  ", or "  ", button to switch


between these two units. Finish by pressing the "  ", button.

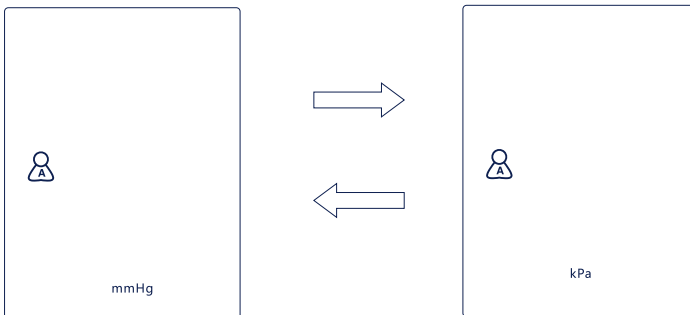
After unit setting, press the "  ", to shut down.

7. USING METHOD OF CUFF


- TIP!: Both arms can be measured.
- Connect the air connector of the cuff to the socket which on the left side of the monitor.

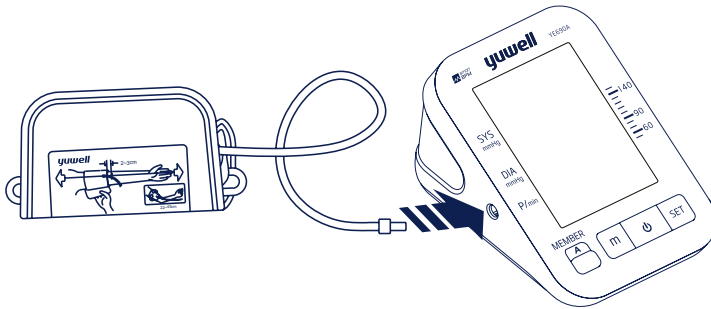
- Place the cuff on table with the stick down. Make the cuff pass through the metal ring, the stick should be showing towards outside.

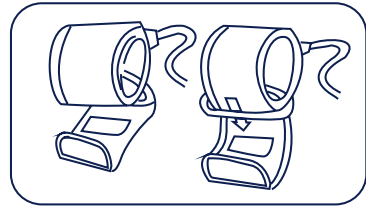- Wind the cuff around the upper arm. (as shown In the picture) Keep the lower edge of the cuff at the position above 2-3cm to the elbow joint and keep the air inlet which insert to the cuff at the inner side of arm.


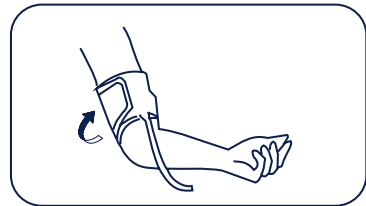
- Stick the cuff tightly.

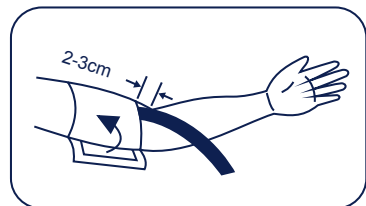
THE CUFF SHOULD BE WEARING COMFORTABLY, AVOIDING TOO TIGHT OR LOOSE.

- Sit straight and flat the arm on table with palm up, also keep the center of cuff and heart at the same level. Also please ensure the air tube not twisted.

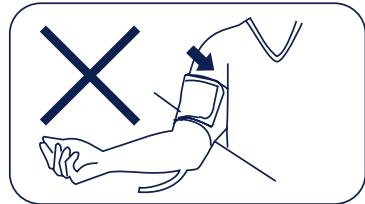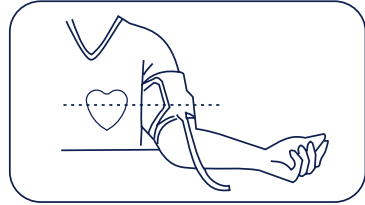8. MEASURING POSTURE

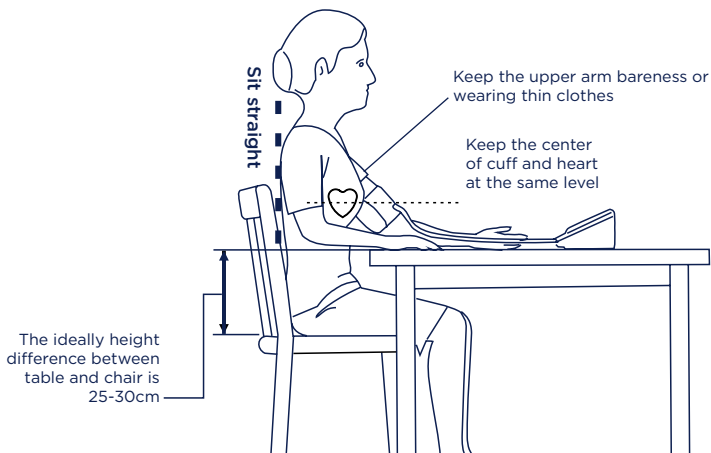8.1. WRONG MEASUREMENT POSTURES

- Do not bend down or body bend forward.
- Do not sit with legs crossed.
- Do not sit on a sofa.(Belly pressure may increase the blood pressure)
- Do not put the arm on the low table.(May increase the blood pressure)

9. START MEASURING

Don't eat, smoke, drink, take bath or do any high impact sport within half an hour before measuring. Measure-


ments shall be taken at the same time every day.

SIT STILL FOR 5 MINUTES BEFORE THE MEASUREMENT

1 Member selecting


Switch the "  ", button to choose Member A or Member B.

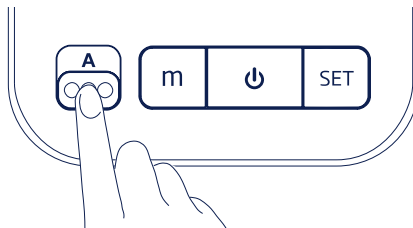
2 Measurement

All icons on screen appear for 1 sec after pressing the "0" button, then disappear. After that the deflation icon flickers which means the system


is in zero testing. Several seconds later the inflation icon flashes which means the zero testing is finished. Then inflation starts.

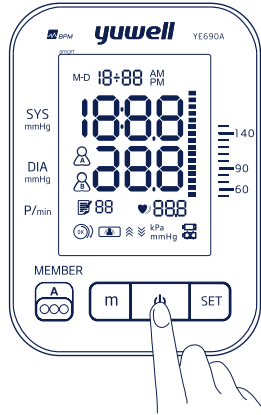Full screen display state

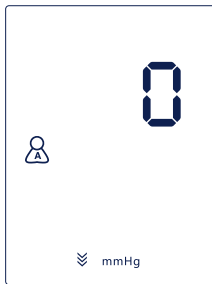The zeroing state

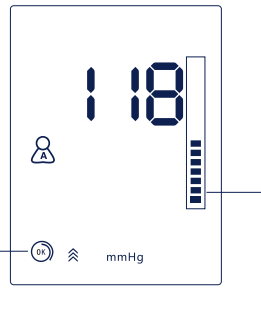
The inflation state

The monitor starts measuring automatically after inflation finished, and the measuring icon „♥“ starts to

flicker, pressure release value starts to decrease gradually.

Please maintain the position during

the measurement and don't speak or move body or hand.

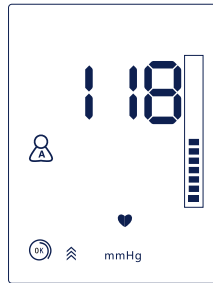The measuring state

3 Measurement finished

After measurement, the monitor will show the SYS, DIA and pulse rate,

then deflate the air automatically.

mmHg display

kPa display


4 Take off the cuff

5 Shut down

Press the „“, will power off automatically in 30 seconds without any operation.

10. BLOOD PRESSURE RANGE INDICATOR

1. If the SYS is higher than 139mmHg or DIA higher than 89mmHg means having high blood pressure. Please contact with doctor for advice.
2. Also the Blood pressure range indicator can show the blood pressure intuitively.

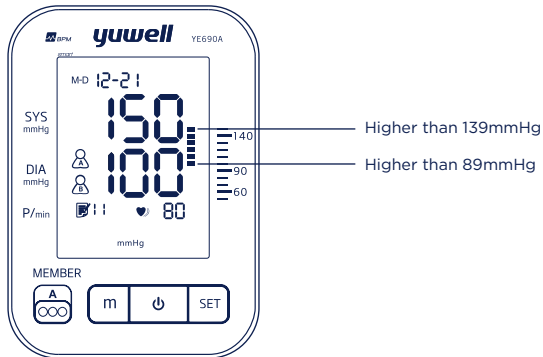- The time interval between two measurements should be at least 2-3 minutes or longer.
- Estimate the blood pressure condition according to the BP classification table, and consult the doctor.

11. MEMORY FUNCTION


The monitor will store the measurement data automatically (including time, date, blood pressure and pul-

se), and the upper limit of records is 74/ user. Press the „m” button to check the records.

11.1. EMERGENCY STOP DURING MEASUREMENT

If the measurement need stopping for uncomfortableness or other reasons, please press the „⏻” button, the measurement will stop immediately

and air releases fast. Take off the cuff manually if the „⏻” button is not working.

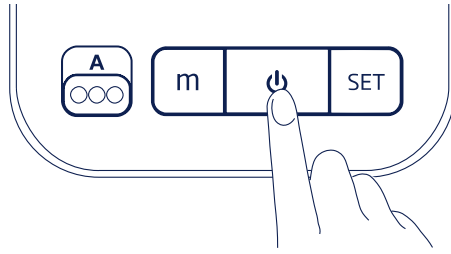
UPPER LIMIT PRESSURE OF AIR INFLATION IS 300MMHG/40.0KPA. DON'T KEEP THE INFLATED STATE FOR A LONG TIME TO AVOID DAMAGE.

- 1 Press the „m“ button for the records.
- Press the ”m„value of the last 3 times measurement.
- B. Reprress the”m„,button to display the 1st group of memory, the serial number is displayed as „1” to „74”. „1” is the latest group and the „74” is the earliest one.

The average value

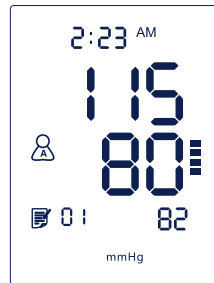The 1st group of memory

The 74th data may be replaced by the 73th when the memory capacity is full.

- Read the recorded data by pressing the „m” button as the sequence: „1”, „2””74”(max). Then return to the 1st. Holding the „m” button to search the
- data quickly.
- Press the „⏻” button to shut off the monitor.

2 Delete the recorded data

Turn the „ „ to member A, then press the both „ „ button and „ „ button until the display shown as the following

picture, which means the recorded data is cleared.

Press the „ „ button to turn off the monitor after deleting.

THIS OPERATION WILL DELETE ALL THE RECORDED DATA.

11. AUXILIARY MEASURING FUNCTION

11.1. CUFF DETECTION

If the cuff winded properly, it will show the „ „ icon. Otherwise it will show the „ „ icon, then press the

„ „ button to stop and wind properly to measure again.


11.2. SWRONG OPERATIO INDICATION

It will show the „ „ icon, if body moves when measuring. Please me-

asure again or it will show inaccurate results.

12. LCD display

- Press the „ „ button and it will show the current time, as shown blew:

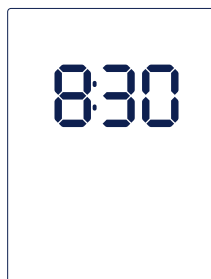- Press the „  „, the current date, time and pressure unit, as shown blew:

Press the „  „, button to shut down the monitor.

12. STATIC MODE

This function is mainly for professional personnel to enter the static

mode to test the monitor through standard pressure gauge.

NORMAL USERS DON'T NEED TO KNOW THIS FUNCTION AND ALSO DO NOT OPERATE. THE COMPANY WILL NOT TAKE ANY RESPONSIBILITY FOR DAMAGE CAUSED BY THIS OPERATION.

12.1. SYSTEM RESTORES

Press the ”  „, button after battery installation, then the screen will show the ”  „, icon, which means the system is in restore testing. Several seconds later, the ”  „, icon disappe-

ars and the air pump starts inflating at the same time, which indicates the test ended. Then press the ”  „, button to stop inflating and take out the batteries to enter the next step.

IT MUST RESTORE THE SYSTEM BEFORE ENTERING THE STATIC MODE, OTHERWISE IT MAY CAUSE INACCURATE RESULTS.


12.2. ENTERING THE STATIC MODE

Press the "  ", button and hold, meanwhile install the batteries. Hold on for about 3 seconds then release the "  ", button. Then screen will

- TIPS: After entering the static mode, if the screen still doesn't show "  ", please operate again as the System restore. Please contact with the local distributor if it still does not work.
- The monitor will automatically power off if there is no operation in 4 minutes.

show the pressure value "  ", the date and time. Now the system has restored and entered the static mode. Now can take the static test.

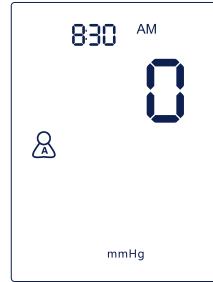The static mode

11. COMMON QUESTIONS OF BLOOD PRESSURE MEASUREMENT

1 What is blood pressure?

Blood pressure is the force exerted by blood against the walls of the arteries. Systolic pressure occurs when the heart contracts. Diastolic pressure occurs when the heart expands. Blood pressure is measured in millimetres of

mercury (mmHg). One's natural blood pressure is represented by the fundamental pressure, which is measured first thing in the morning while one is still at rest and before eating.

2 What is Hypertension and how is it controlled?

Hypertension, an abnormally high arterial blood pressure, if left unattended can cause many health problems including stroke and heart attack. Hypertension can be controlled by altering lifestyle, avoiding stress, and with

medication under a doctor's supervision.

To prevent Hypertension or keep it under control: Do not smoke, exercise regularly, reduce salt and fat intake, have regular physical checkups, maintain proper weight.

3 Why measure Blood Pressure at home?

Blood pressure measured at a clinic or doctor's office may cause apprehension and can produce an elevated reading, 25-30 mmHg higher than that measured at home. Home me-

asurement reduces the effects of outside influences on blood pressure readings, supplements the doctor's readings and provides a more accurate, complete blood pressure history.

4 WHO Blood Pressure Classification?

Standards to assess high blood pressure, without regard to age, have been established by

the world Health Organization (WHO), as shown below:

RANGE	SYSTOLIC PRESSURE KPA/MMHG	DIASTOLIC PRESSURE KPA/MMHG	COUNTER-MEASURES
Hypotension	<12.0/90	≤8,0/60	Consult dr.
Ortho-arteriotony	12.0/90~18.5/139	8.0/60~11.9/89	Self check
Mild hypertension	18.7/140~21.2/159	12,0/90~13,3/100	Consult dr.
Medium hypertension	21.3/160~23.9/179	13.3/100~14.5/109	Consult dr.
Severe hypertension	≥24.0/180	≥14,7/110	Danger! Go to hospital as soon as possible

THERE IS NO DEFINITION ABOUT HYPOPIESIA, AND GENERALLY SYS (SYSTOLIC PRESSURE) LESS THAN 90MMHG OR DIA (DIASTOLIC PRESSURE) LESS THAN 60MMHG IS CALLED HYPOTENSION.

5 Blood pressure variations?

An individual's blood pressure varies greatly on a daily and seasonal basis. It may vary by 30 to 50 mmHg due to various conditions during the day. In hypertensive individuals variations are even more pronounced. Normally, the blood pressure rises while

at work or play and falls to its lowest levels during sleep. So do not be overly concerned by the results of one measurement. Take measurements at the same time every day using the procedure described in this manual to get to know your normal blood

pressure. Regular readings give a more comprehensive blood pressure history. Be sure to note date and time when recording

your blood pressure. Consult doctor to interpret the blood pressure data.

6 When is the best time to measure blood pressure?

- After urination, before breakfast in the morning.
- Before sleeping at night.
- Before taking medicine.
- Please keep a stable body state and mind every time measuring. We suggest taking measurements at a same time every day.

12. COMMON FAULT AND TROUBLE SHOOTING

COMMON FAULT	SOLUTIONS
It doesn't work after pressing the ' ' button with batteries installation	Check the batteries are installed correctly
	Replace new batteries
Multiple occurrences of measuring failure, or measured value is low (or high)	Check the connection and winding of cuff
	Check if the cuff winded too tight or too loose. Take off your clothes if rolled too tight
	Please ensure a quiet, relaxed body state o. Deep breathing to relax yourself before measurement
The monitor is in good condition, but the each measuring result are different	Please read the „blood pressure variations“ carefully
The value is different from that measured at a clinic or doctor's office	Write down the value every day, and consult a doctor
Pump works, but the pressure doesn't rise	Check whether the cuff has connected well

The table below shows the possible fault displaying icon, possible reason, and solutions.

WRONG INDICATION	FAULT CAUSE	SOLUTION
Err4	Unable to measure pressure	Fasten cuff correctly before measurement
Err5	Pressurizing error	Check if there is air leakage from the cuff
Err6	Pressurizing error caused by arm or body motion	Keep arm and body still and measure again
Err7	Cuff is too loose or fall off	Fasten cuff tightly
Err8	Pressure exceeds the maximum value (300mmHg)	Measure again please
Showing low voltage icon „  „	Battery is low	Replace new batterie
howing the „  „icon	Battery is running out	Replace new batteries

IF THE SITUATIONS CANNOT BE SOLVED OR UNEXPECTED PROBLEM HAPPENS, PLEASE CONSULT THE LOCAL DISTRIBUTOR.

13. MAINTENANCE

Please observe the following items to protect the device and ensure the accuracy of measurement.

PLEASE STORE THE MONITOR AND ACCESSORIES PROPERLY AFTER USE.

DO NOT PLACE THE MONITOR AND ACCESSORIES IN HIGH TEMPERATURE, MOISTURE, DUST, OR EXPOSURE TO SUNSHINE.

THE CUFF CONTAINS AN AIRBAG INSIDE, PLEASE CARE IN APPLICATIONS, DO NOT FOLD, PULL OR TWIST IT.

DO NOT DISASSEMBLE OR REPAIR THE DEVICE WITHOUT AUTHORIZATION.

DO NOT REPLACE THE PARTS WITHOUT AUTHORIZATION.

PLEASE CLEAN THE MONITOR WITH SOFT DRY CLOTH. IF IT'S NECESSARY, PLEASE USE WIPED SOFT CLOTH WITH WATER OR NEUTRAL DETERGENT BEFORE CLEANING BY SOFT DRY CLOTH.

USING ABSORBENT COTTON TO WIPE GENTLY WITH RUBBING ALCOHOL, DISINFECTION OF THE MACHINE WHEN IT NECESSARY. DO NOT USE DETERGENT TO CLEAN.

DO NOT LET WATER SEEP INTO THE DEVICE.

TIPS: We advice to calibrate the monitor according to local laws and regulations.

13.1. FEATURES AND TECHNICAL PARAMETERS

1. Features:
 - Large-scale backlight LCD
 - Cuff detection
 - IHB(Irregular heartbeat) indicator
 - Misoperation detection
 - 2X74 groups of memory
 - Average value display of measurement
 - Intelligent pressurized
 - Time and date display
 2. Technical parameters
 - Displaying: LCD digital display
 - Measuring method: oscillation mensuration
 - Measuring range:
 - Systolic pressure: 60-230mmHg (8.0-30.7kPa)
 - Diastolic pressure: 40-130mmHg (5.3-17.3kPa)
 - Pulse rate: 40-200 time/min
 - Precision:
 - Pressure: within 3mm g(0.4kPa)
 - Pulse rate: within 5% of reading value
 - Working system: continuous operation
 - Electric classification: Class II, type BF application part
 - Pressure sensor: Semiconductor pressure sensor
 - Pressurization: Automatic pressurize Depressurization: Automatic air releasing Power supply: 4 AA batteries or 6V/ 600mA DC power
 - Maximum battery life: 300 times
 - Suitable upper-arm circum-
- ference: 22-45cm Dimension: 133x100x83(mm)
- Weight: 351.8g
 - Operation and storage conditions:
 - Working temperature: +5°C-+40°C
 - Relative humidity: 15%-80%
 - Working atmosphere pressure: 80kPa-105kPa
 - Transport and storage temperature: -20°C~ +55°C
 - Transport and storage relative humidity:15%-80%, no condensation
 - Transportation and storage atmosphere pressure: 80kPa ~ 105kPa
 - Operation environment: Avoid electromagnetic interference, violent shock and noise environment. Radio frequency interference

13.2. FEATURES AND TECHNICAL PARAMETERS

There is the potential risk of radio frequency interference between the device and other devices. If there is, please find out the problems and take the following measures:

1. Turn off the device, and turn on again,
2. Change the direction of the device,
3. Keep the product away from the interferential devices.

smart
2M BPM

JIANGSU YUYUE MEDICAL
EQUIPMENT & SUPPLY CO., LTD.
Yunyang Industrial Park
2120300 DanYang, Jiangsu
PEOPLE'S REPUBLIC OF CHINA

Shanghai International Holding
Corporation GmbH
Eiffestraße 80
20537 Hamburg, Niemcy

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